



Ensure Patient Safety  
Protect your Investment

## DP360

Premium Steam Distillation Appliance

### The Best Steam-Distilled Water Solution for Dentists

**PROBLEM #1: BACTERIA & BIOFILM.** According to the American Dental Association, dental unit waterlines have been shown to harbor, in significant numbers, a wide variety of bacteria, fungi, and protozoans.

**PROBLEM #2: VERY SENSITIVE AUTOCLAVES.** New autoclaves and Statim machines require pure distilled water on a consistent basis or they will not operate properly and their lifespan could be shortened.

**SOLUTION:** Dentists should use only pure, steam-distilled water in their Dental Units and autoclaves.

### Recommended by Unthank Design Group

*"For patient safety, we have long advocated using steam-distilled water in a closed water system for dental delivery. A DENTAL PURITY in-office distillation system provides for a convenient continuous supply of the purest water."*



Michael Unthank DDS  
Registered Professional Architect & Dentist  
Unthank Design Group ([www.unthank.com](http://www.unthank.com))

## DP360

### Premium Steam Distillation Appliance

Our DP360 Steam-Distillation Appliance is perfect for smaller dental offices. Simply fill up the boiling tank with tap water, turn the appliance on and the DP360 will produce about 3 liters of distilled water per cycle.

- Manual operation (operates much like a coffee maker).
- No installation required.
- Simple maintenance.
- High-quality, stainless steel construction.
- Perfect for small to mid-sized dental offices.



Visit us on-line at [www.DentalPurity.com](http://www.DentalPurity.com)

Or call one of our customer care representatives today at **1-800-617-2401.**

Dental Purity is a Brand of Pure & Secure, LLC. 4120 NW 44th Street, Lincoln, NE 68524,

## PROVIDE SAFE, PURE WATER FOR YOUR PATIENTS

Steam-distilled water is an important key in preventing bacteria and biofilm build-up in your water lines so you can provide a safe, pure water supply for your patients.

## TAKE CONTROL OF YOUR DISTILLED WATER

The DP360 allows you to take control of your own supply of pure, steam-distilled water. No more running out of distilled water, last minute trips to the store or the wasted space of bottled water.

## SIMPLE OPERATION

The DP360 is as simple to operate as a coffee maker. Simply fill the boiling tank with tap water, turn the unit on and it will produce 3 liters of distilled water in about three hours.

## NO INSTALLATION REQUIRED

The DP360 does not require installation. Just plug it in and it's ready to go.

## 40 YEARS OF TRUST

For 40 years, dentists around the world have used our premium steam distillation appliances to ensure that they are using the highest quality water at all times.

## FAR SUPERIOR TO COMPETITIVE DISTILLERS

The DP360 is far superior to competitive distillers on the market. Here's why...

**QUALITY CONSTRUCTION** The DP360 is constructed of stainless steel and top-quality components to provide years of dependable service.

**LONG LIFESPAN** The DP360 is a heavy duty appliance and is built-to-last. The DP360 has a proven track record and a 15-year limited warranty.


**SIMPLE MAINTENANCE** Since the DP360 never boils dry it is much easier to maintain than competitive distillers. Simply dump the remaining water out of the boiling tank each time you use it in order to greatly reduce the amount of scale build-up in the boiling tank.

**AN EXCELLENT INVESTMENT** When comparing the DP360 to other manual distillers, which have a very short lifespan, or to bottled water it's clear that the DP360 is an excellent investment that will pay you back many times over.

## SPECIFICATIONS

item #.....	DP360
width.....	13.75" (35 cm)
depth.....	9.875" (25 cm)
height (120V).....	15.75" (40 cm)
height (240V).....	16.5" (42 cm)
shipping weight.....	32 lbs. (14.5 kg)
machine weight.....	26 lbs. (12 kg)
current draw (120V).....	6.9 amps
current draw (240V).....	3.5 amps
total wattage.....	800 watts

Specifications and product list may be subject to change without notice.

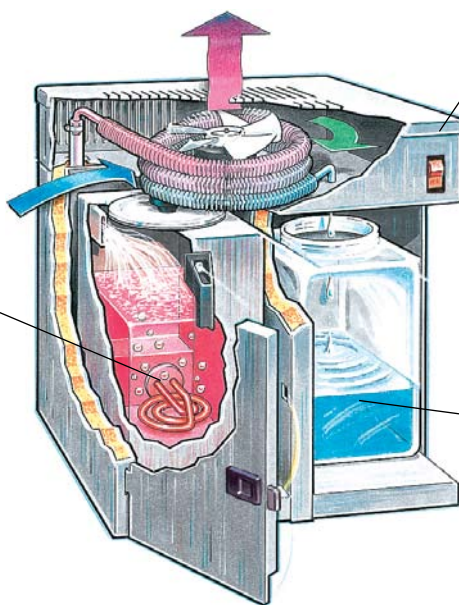

**TIMER SHUT-OFF** The DP360 is equipped with a unique timer shut-off switch that stops the DP360 at the end of each distilling cycle so that it doesn't boil dry. Other manual distillers on the other hand, turn off by boiling dry and then superheating, which causes an overheat switch to turn the machine off. This method of turning off can create a fire hazard (there have been large recalls of competitive distillers for this very reason).

**CONSISTENT PURITY** When other distillers boil dry the temperature of the heating element superheats which can cause unwanted contaminants to vaporize and contaminate the distilled water. The DP360 does not boil dry so it doesn't have this problem.